
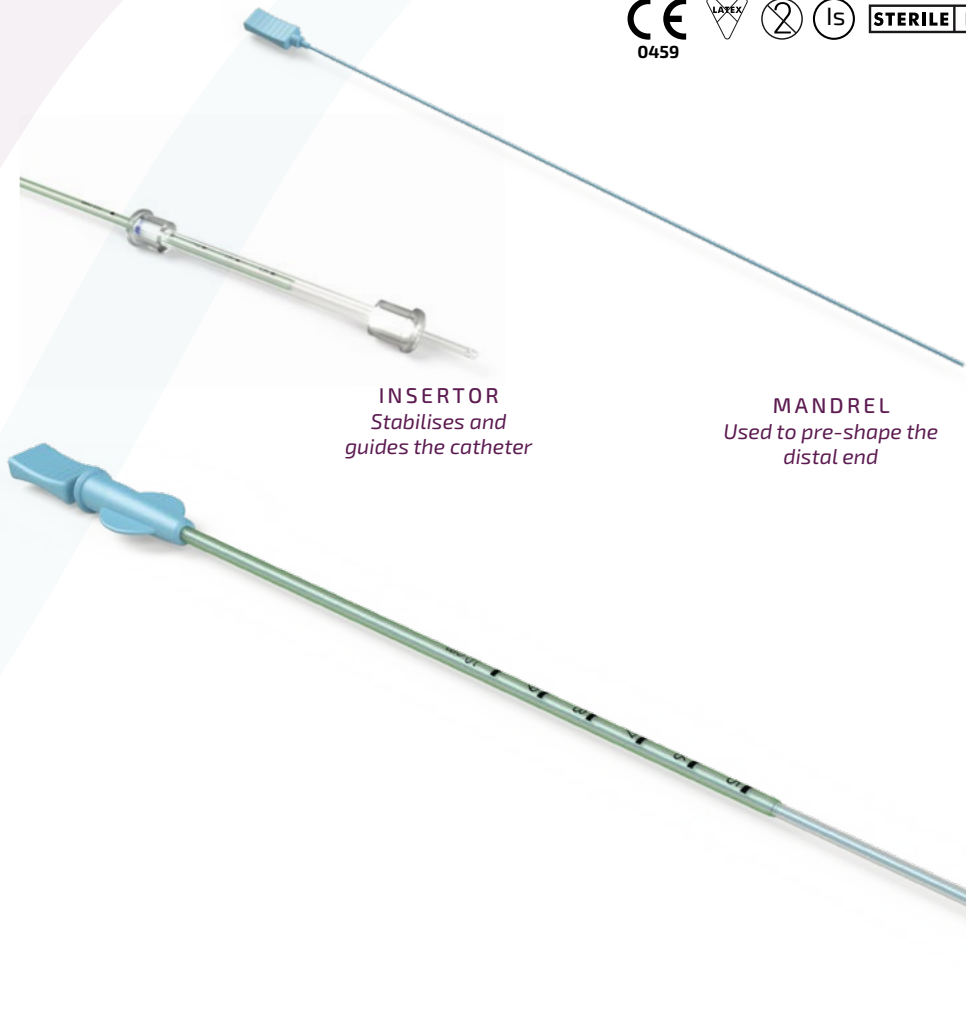


## PM-LIFE® IU

## INSEMINATION CATHETER

31

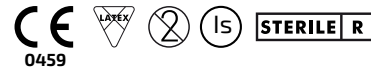

**INSERTOR**  
Stabilises and  
guides the catheter

**MANDREL**  
Used to pre-shape the  
distal end

The PM-LIFE® IU intrauterine insemination catheter is used to perform insemination during medically assisted reproduction (MAR), by introducing washed and prepared sperm into the uterine cavity. Some product references feature an insertor, to stabilise and guide the catheter during insertion, and/or a mandrel to pre-shape the distal end of the catheter in order to adapt to the patient's morphology.

**ADVANTAGES**

- Centimetric graduations for better control of catheter insertion.

**AVAILABLE REFERENCE(S)**

Reference	Diameter for catheter (mm / Fr)	Inner diameter volume (mm / ml)	2 lateral holes diameter (mm)	Useful length (mm)	Length of the distal part (mm)	Insertor	Mandrel
PS1260500	1.65 / 5	1.1 / 0.20	1	175	45	-	-
PS1270500	1.65 / 5	1.1 / 0.20	1	175	45	Yes	-
PS4130500	1.65 / 5	1.1 / 0.20	1	175	45	-	Yes
PS4160500	1.65 / 5	1.1 / 0.20	1	175	45	Yes	Yes
PS1260600	2 / 6	1.6 / 0.40	1.3	175	45	-	-
PS1270600	2 / 6	1.6 / 0.40	1.3	175	45	Yes	-
PS4130600	2 / 6	1.6 / 0.40	1.3	175	45	-	Yes
PS4160600	2 / 6	1.6 / 0.40	1.3	175	45	Yes	Yes

**ABOUT THIS MEDICAL DEVICE**

1 <sup>st</sup> release date	2002
Single use	Yes
Class	Is
Sterile	Yes
Sterilisation method	Sterilised via gamma radiation
Contains latex	No
Lifetime	5 years
Certification	CE
Packaging	Box of 25 units in individual sterile wrapping
Storage	In its original box in a cool and dry place

PRINCE MEDICAL SAS  
ZA La Sente Du Moulin  
64 rue Gutenberg - 60530 Ercuis - France  
www.prince-medical.com  
Tel: +33 (0)3 44 26 54 76  
prince-medical@omerin.com

The information provided on this data sheet is indicative and subject to change without notice. PRINCE MEDICAL SAS shall not be held responsible for any incidents resulting from inappropriate use, failure to comply with the rules and regulations in force, or failure to follow the instructions for use provided. Medical devices sold by PRINCE MEDICAL SAS should only be used by a duly qualified practitioner or authorised person.

® Registered trademark of PRINCE MEDICAL SAS. Drawings and photos are not contractually binding. Reproduction is prohibited without the prior consent of PRINCE MEDICAL SAS.

PRINCE MEDICAL  
DEVICES